




G8 HPLC Reporting Software

Reporting Software Expands Functionality for Tosoh G8 HPLC Analyzer

- User Friendly
- Data Management
- Levey-Jennings
- Chartable Reports
- Lot # Documentation

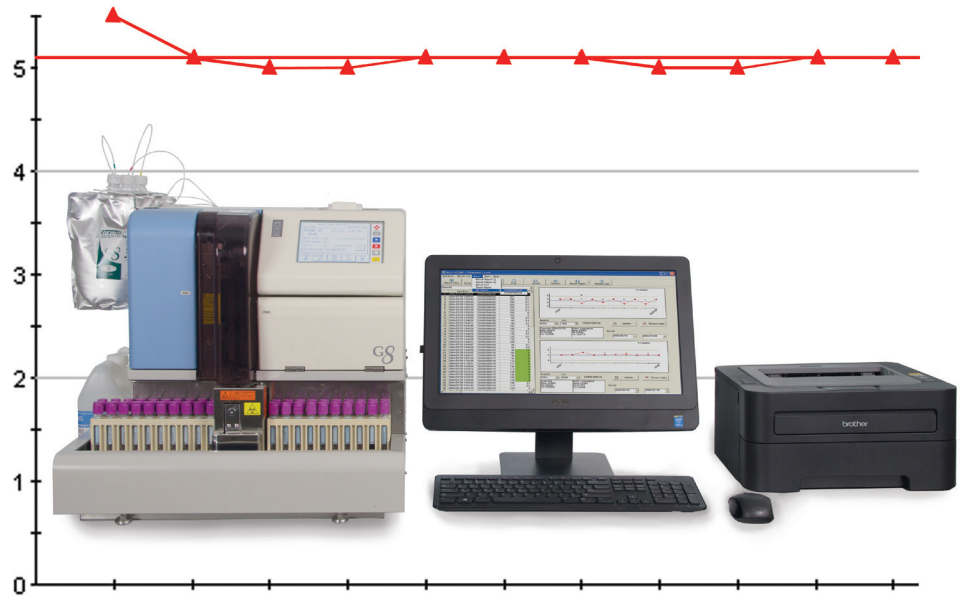

TOSOH AUTOMATED GLYCOHEMOGLOBIN ANALYZER HLC-723G8

G8


TOSOH BIOSCIENCE

G8 HPLC Analyzer Features


- G8 HPLC Analyzer with touch screen - Simplicity of operation
- Cation exchange chromatography - Gold Standard method used for DCCT
- SA1c in 1.6 minutes - Quick turn around time to optimize patient care
- Multifunctional STAT sample position - No dilutions - Saves time
- Physical design - Compact size fits anywhere
- Mechanical design - Low CVs for a high level of precision
- NGSP certified - National Standardization
- IFCC certified - International Certification

G8 Reporting Software Features


- Customizable G8 database management - Flexibility in result reports
- Automatic printing of patient demographics with chromatogram - Provides chartable reports
- Trend analysis charts - Easy review of patient SA1c results over time
- Result tracking by lot - Enhanced documentation for ease of inspections

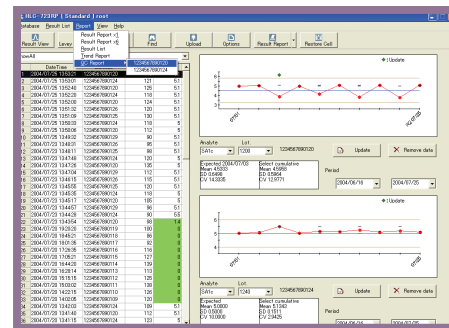


Levey-Jennings

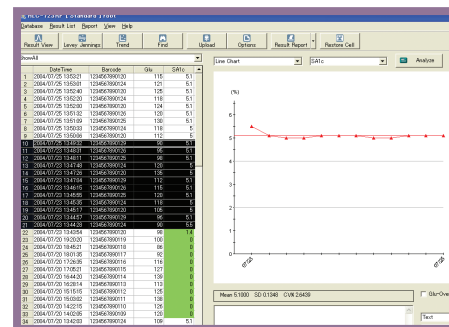

Trend

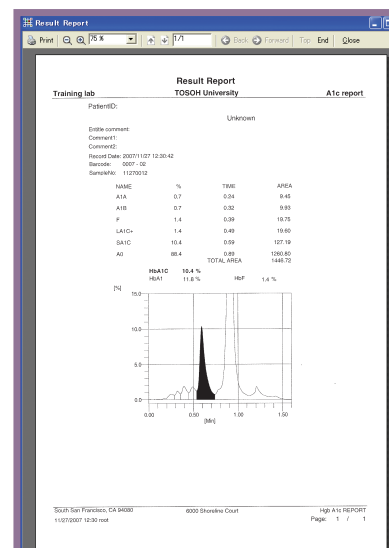

Result Report - 1 Chromatogram

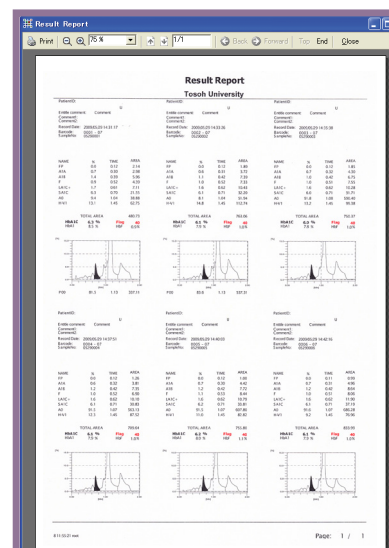

Result Report - 6 Chromatogram



G8 HPLC Reporting System

Complete HPLC Reporting System Includes:

- | | |
|---|--|
| G8 HPLC Analyzer (Instrument) | Printer |
| G8 Reporting Software | RS 232 Cable 9PINF to 9PINF |
| Dedicated PC | Serial Ports/Serial Port Cables |
| OptiPlex 3030 All-In-One Desktop | 1 - If Stand Alone |
| Windows 7 Professional (Operating System) | 2 - If the PC is interfaced to Host Computer |
| Intel® Core™ i3-4130 Processor (3.4 GHz) | Wireless Keyboard |
| Memory - 4 GB | Wireless Mouse |
| Hard Drive - 500 GB | |

For more information call (800) 248-6764
or visit www.tosohbioscience.us



Tosoh Bioscience, Inc.
6000 Shoreline Court, Suite 101
South San Francisco, CA 94080
Tel: (800) 248-6764
Fax: (650) 615-0415
www.tosohbioscience.us